


# DIAGNOSTIC AND LABORATORY EQUIPMENT



**T**

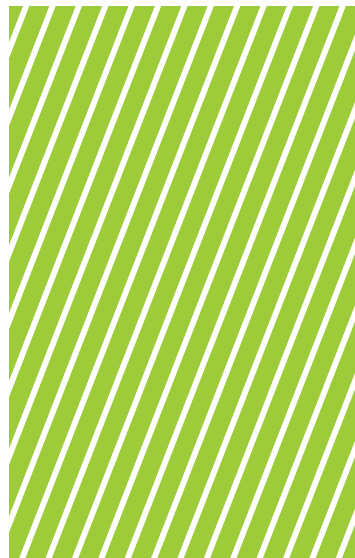

**DEMONSTRATION MODELS FOR FIXED ORTHODONTICS**

Two-tone acrylic normo-occlusion demonstration models are useful for the visualization of the bracket characteristics for both patients and potential customers. Available with brackets of various types, steel, ceramics and even self-ligating; single Low Profile tubes are present on the molars in all demonstration models.

Pack of 1



**T3472-00**  
**DEMONSTRATION MODEL - MINI DIAGONALI®**



**T3442-00**  
**DEMONSTRATION MODEL - MIDI DIAGONALI**



**T3461-00**  
**DEMONSTRATION MODEL - STEP SYSTEM 2.0**





**T3457-00**  
**DEMONSTRATION MODEL - AQUA SL**




**T3411-00**  
**DEMONSTRATION MODEL - INTERACTIVE SL**



**T3456-00**  
**DEMONSTRATION MODEL - GHIACCIO**



**T3458-00**  
**DEMONSTRATION MODEL - AQUA**

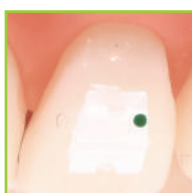

**T3061-00**  
**PLASTIC**  
**MODEL BASES**

Impact resistant, white plastic model bases, with rear hinges to maintain the occlusion.

Square-edged according to international standards and a constant height that permits to make models comparable over time and in all views for the correct diagnosis and for monitoring during therapy. For use with the practical positioning plate T3062-02 supplied with the package.

Precision and ease of working models, that can also be performed directly in the dental office.

**Pack content:**

- 100 upper bases
- 100 lower bases
- 2 position plates with guide rails and labels



**T3062-00**  
**SPARE HINGES**  
**FOR MODEL BASES**

Pack of 10

**T3062-01**  
**SPARE LABELS**  
**FOR MODEL BASES**

Pack of 100

**T3062-02**  
**POSITION PLATES**  
**FOR MODEL BASES**

Pack of 10 with 20 guide-rails

**T1121-00**  
**VERTICAL MODEL**  
**OCCLUDER**

Entirely made of stainless steel with high precision. The double sliding guides prevent from any undesired lateral movement. Models can be disengaged from the occluder by loosening the locking screw.

**Locking screw T1122-00.**





**RUBBER  
MODEL BASE  
FORMERS**

**T3320-00** Large  
Pack of 2,  
1 upper and 1 lower

**T3321-00** Small  
Pack of 2,  
1 upper and 1 lower



**T3355-00  
MODELBOX**

Transparent PVC model containers supplied with pre shaped sponges for full protection of models.

**Pack of 100**



**T3529-00  
ACRYLIC TEETH  
FOR TYPODONT T3330-00**

Purposely designed for Typodont T3330-00, they are necessary for making practice with the direct bonding brackets.

**Pack of 28 teeth**



**WAX BASES FOR TYPODONT**

- T3332-00** Class I
- T3336-00** Normal occlusion
- T3337-00** Class II, I div.
- T3338-00** Class II, II div.
- T3339-00** Class III

**Pack of 2,**  
upper and lower





**T1315-00**  
**LEONE HYDROFLASK**

Orthodontic cold curing acrylics require the polymerization under pressure technique. This kind of polymerization enables the vaporization of monomer and condensing of acrylic, which will perfectly fit the plaster model. Plates are coming out perfectly finished and can be immediately polished. The Leone hydroflask allows you to make use of the polymerization under pressure technique easily and economically.

**Technical features:**

Inside dimensions: 100 x 115 mm  
Outside dimensions: 111 x 130 mm  
Pressure gauge  
Exhaust valve  
Made of hard aluminium  
Built-in thermometer



**T1310-00**  
**LEONE HYDROFLASK**

With the same features as T1315-00, without thermometer and slightly smaller.

**Technical features:**

Inside dimensions: 100 x 88 mm  
Outside dimensions: 111 x 105 mm  
Pressure gauge  
Exhaust valve  
Made of hard aluminium

**REPLACEMENT PARTS  
FOR HYDROFLASK**

<b>T2300-00</b>	Hydroflask lid
<b>T2301-00</b>	Rubber seal for lid
<b>T2302-00</b>	Air valve
<b>T2303-00</b>	Thermometer for T1315-00
<b>T2304-00</b>	Pressure gauge



**T1311-00**  
**PRESS FOR HYDROFLASK**

Convenient press for hydroflasks T1310-00 and T1315-00, made of pressure cast aluminium. Extremely easy to use.

**TRAYS  
FOR MODEL BASES**

<b>T3360-00</b>	Blue
<b>T3361-00</b>	Red

Two drawers for storing six model bases each box.

Outer sizes: 250 x 200 x 80 mm.

**Pack of 16 boxes**

**TRAYS  
FOR STUDY MODELS**

<b>T3350-00</b>	Blue
<b>T3351-00</b>	Red

Two drawers for storing eight models each box.

Outer sizes: 285 x 190 x 55 mm.

**Pack of 24 boxes**





**T3603-00**  
**ACETATE SHEETS FOR**  
**CEPHALOMETRIC TRACING**

Sheets with double matt and polish surface. Sizes: 203 x 254 mm.  
**Pack of 100 sheets**



**T3608-00**  
**MARKING**  
**PENCIL WHITE**

For marking arch wires, bite planes, etc.  
**Pack of 1**



**T3612-01**  
**CEPHALOMETRIC**  
**PROTRACTOR**

The angular measurement is facilitated by the ruler; the cut-out tooth symbols of incisors, molars and circular holes assist in the anatomical tracing over the x-ray film.  
**Pack of 1**



**T1565-00**  
**SAW**



**BLADES FOR SAW**  
**T1565-00**

**T1566-00** 12 lame per gesso

**T1567-00** 12 lame per metallo



**SAW**  
**FOR SET-UP**

**T1561-00** Saw with 1 blade

**T1563-00** 12 large serrated blades

Easy change of blades and perfect handiness.

**LEONE COPOLYMER DUPLICATION PLATES**

Sizes: 130 x 140 mm, 2 or 3 mm thick, these plates are made of a special material getting soft and gelatinous at a temperature of 90° C. Just dip in boiling water, lay on the model and press over with the duplicating press: a fine detailed model impression is obtained.

**T0900-02** Pack of 50 plates 2 mm thick

**T0900-03** Pack of 50 plates 3 mm thick



**T0900-92** Pack of 500 plates 2 mm thick

**T0900-93** Pack of 500 plates 3 mm thick



**T1200-00**  
**PLIER RACK**



**T1203-00**  
**PLIER RACK**



**T1212-00**  
**STANDARD ARCHWIRE HOLDER**